


# **PRONTO<sup>®</sup> V<sub>3</sub>**

*Extraction Catheter*

Specially designed for percutaneous thrombus aspiration

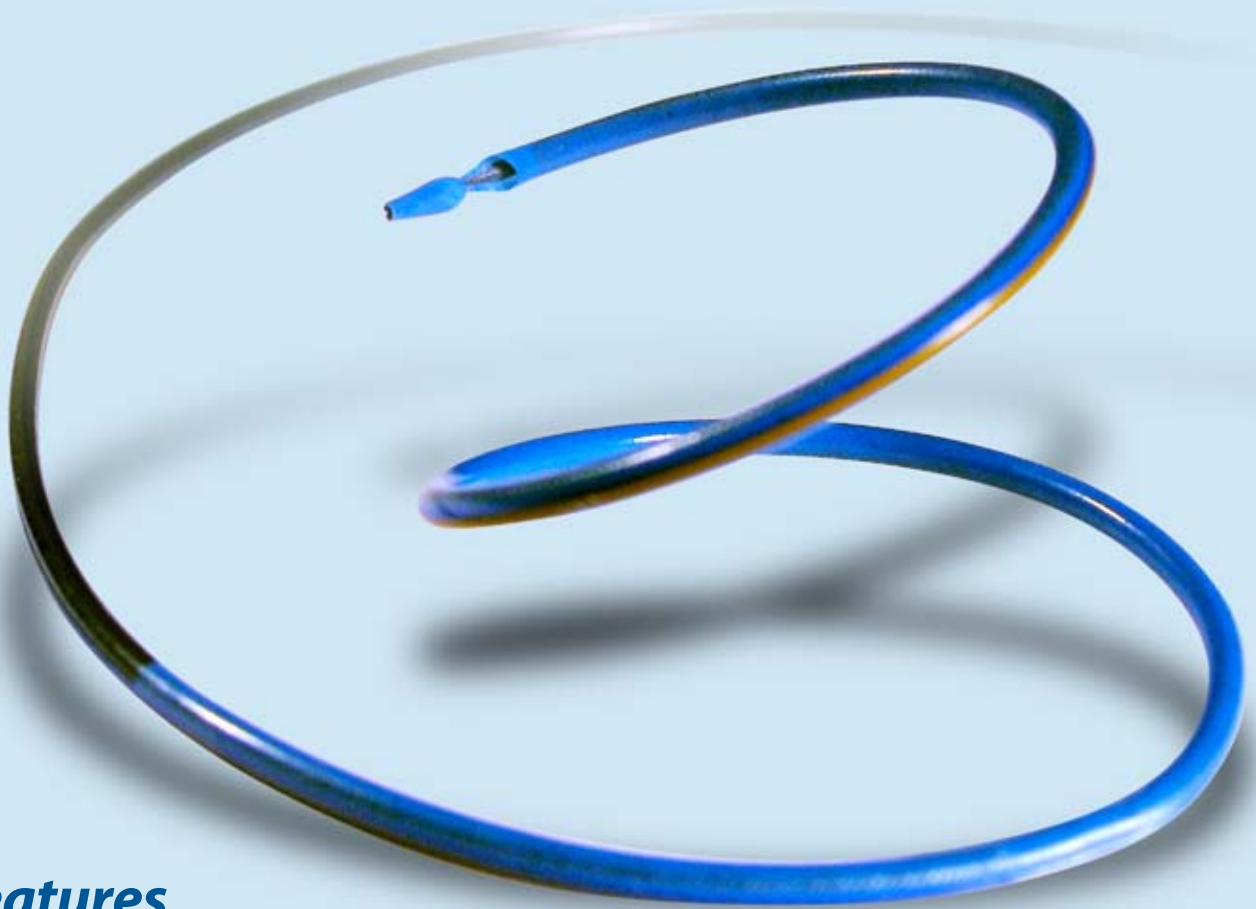

## **V<sub>3</sub> Features**

- Full-length wire braid
- Uniformly large extraction lumen
- Hydrophilic coating
- Lubricious wire lumen
- Tapered distal tip

# The Next Generation of Soft Thrombus Extraction ...

## **PRONTO!**

*Specially designed  
for percutaneous  
thrombus aspiration*



### *Simple one-operator deployment*

- Rapid exchange 140cm catheter
- Positioning marks at 95cm (single) and 105cm (double)
- Two locking 30ml aspiration syringes
- 6F ( $\geq 0.070$ " I.D.) guide compatible

The **Pronto V3** extraction catheter is indicated for:

- the removal of fresh, soft emboli and thrombi from vessels in the coronary and peripheral vasculature

The **Pronto V3** extraction catheter is contraindicated in:

- vessels <2mm in diameter
- the removal of fibrous, adherent or calcified material (e.g., chronic clot, atherosclerotic plaque)
- the venous system
- cerebral vasculature

V3


# The V<sub>3</sub> Advantages ...

## Advanced Catheter Design

- Fully braided shaft increases flexibility, pushability and kink resistance
- Hydrophilic coating for smooth navigation
- Single extrusion for uniformly large extraction lumen diameter
- Lubricious wire lumen




## Uniformly Large Extraction Lumen

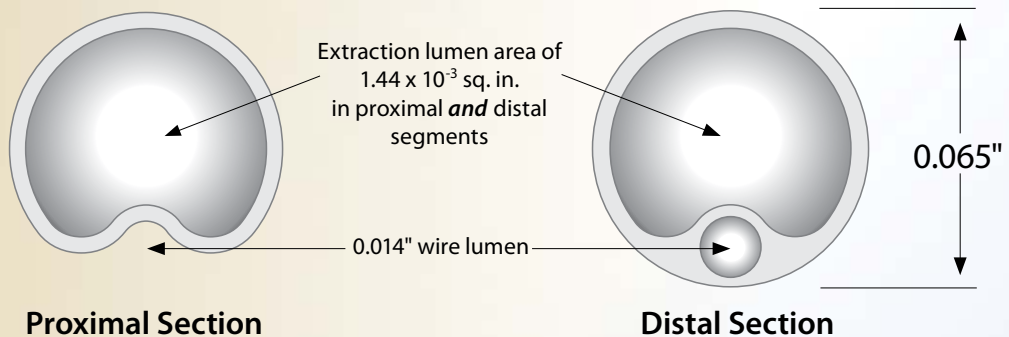

## Tapered Silva Tip\*

- Tapered entry profile—*designed to improve navigation*
- Self-centering—*designed to extract the thrombus while preventing vessel wall adhesion*
- Sloped extraction lumen opening—*designed to accommodate larger thrombi*

### Pronto V<sub>3</sub> Silva tip

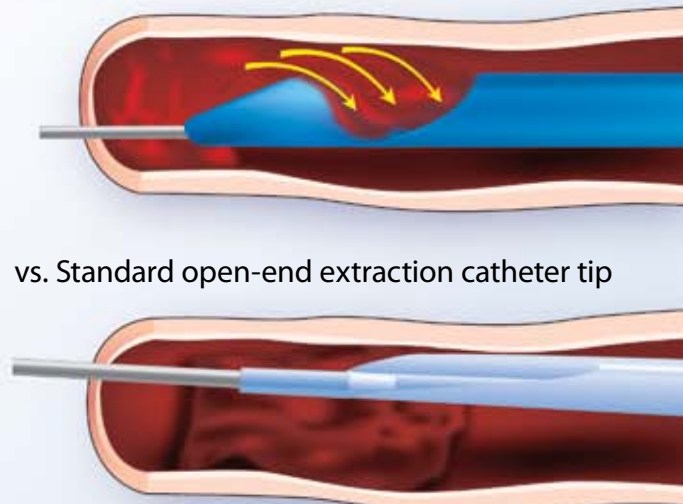

# PRONTO® V3

## Extraction Catheter

The Pronto V3 extraction catheter is indicated for the removal of fresh, soft emboli and thrombi from vessels in the coronary and peripheral vasculature

Model	Description	Box Quantity
5003	Pronto V3 Extraction Catheter	1



### Specifications:

Guide catheter compatibility:	≥6F (0.070"/1.78mm min. I.D.)
Guidewire compatibility:	≤0.014"
Working length:	140cm
Positioning marks:	95cm = single white mark 105cm = double white marks
Rapid exchange section:	20cm
Marker band:	Located ≈2mm from the distal tip
Syringes:	Two locking 30ml polycarbonate
Filter basket:	70µm mesh

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

U.S. and International patents pending.  
Pronto is a registered trademark of Vascular Solutions, Inc.

©2012 Vascular Solutions, Inc. All rights reserved. ML1559 Rev. J 04/12



**Vascular Solutions, Inc.**  
6464 Sycamore Court  
Minneapolis, Minnesota 55369 USA

Customer Service:  
United States: 888.240.6001  
International: (001) 763.656.4298  
customerservice@vasc.com

[www.vasc.com](http://www.vasc.com)